

Invacare® | **Action<sup>4</sup><sub>NG</sub>**



COMPATIBLE



**Where stability meets style!**

A comfortable, foldable, lightweight wheelchair, the *Invacare Action<sup>4</sup><sub>NG</sub>* is ideally suited for everyday use.

A key feature of this next-generation model is the innovative Invacare hallmark of cross-compatibility, a shared platform with other models in the *Action<sub>NG</sub>* Family as well as other wheelchairs in our range means that users and advisors can easily update or refurbish their wheelchairs and enjoy maximum opportunities for modularity and personalized adjustment.



**Yes, you can.®**



### Easy rolling and maneuverability

The particularities of the *Invacare Action<sup>4</sup><sub>NG</sub>* frame lie in its stability, resulting in easier wheeling and better utilization of the user's body movements. With its main balance options it is easy to adapt the *Action<sup>4</sup><sub>NG</sub>* to other types of wheeling characteristics. The seat has optional adaptable seat depth and seat angle to ensure that each user obtains the necessary postural stability and support. The *Action<sup>4</sup><sub>NG</sub>* has the characteristic *Action<sub>NG</sub>* design bend in the lower frame.

### Compatibility with other Invacare wheelchairs

Almost all parts of the *Action<sup>4</sup><sub>NG</sub>* frame are interchangeable with other members of the *Action<sub>NG</sub>* Family and even with other Invacare wheelchairs, including *Rea<sup>®</sup>*. All parts such as armrests, legrests, wheels etc, are compatible with each other.

### Easy transport with low weight

By removing the rear wheels and folding the frame, the *Action<sup>4</sup><sub>NG</sub>* easily slides in behind the front seat or slips into the boot of a car. Without rear wheels,

the *Action<sup>4</sup><sub>NG</sub>* weighs 8.3 kg. Simple handling and folding makes for great ease of transport.

## Features and Options



#### Adjustable backrest upholstery

As an option, the backrest can be equipped with tension-adjustable straps, shaping the backrest to individual user needs.



#### Seat-depth adjustability

Decrease seat depth by 50 mm, to optimize seat size.



#### Backrest bar

With tension-adjustable upholstery a backrest bar is required for extra stability.



#### Drip stand

## Features and Options




### Legrests assortment

① Elevating legrest from 0° to 80°, ② fixed 80° legrest and ③ amputee support are available with either ④ fixed, ⑤ angle and depth adjustable or ⑥ one piece footplate.



### Optional rear wheel positions

Active position plate offers better balancing possibilities for superior driving performance.

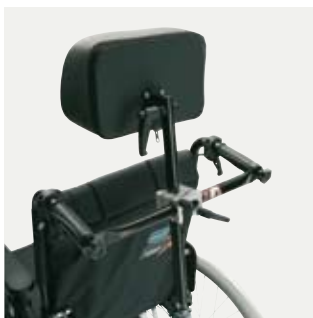

### Headrest

Can be adjusted on the backbrace in width and depth to provide optimum positioning.



### Hub brakes

Hub brakes available as an option, helping assistants to stop the wheelchair and reduce speed.



### Anti-tip device

Prevents wheelchair from tipping over and is height adjustable to suit the different seat heights available.



### Amputee support

Height, angle and depth adjustable, the amputee support is also swing away.



### Heavy Duty Version up to 160 kg

A wider choice of seat widths (505, 555 and 605 mm), reinforced chassis and a carrying capacity of 160 kg.



### Hemi version

A 385 mm seat height, dual handrim propulsion and hemi armrest are available for hemi users.



**Technical data**

<b>Self propulsion</b>	380 / 405 / 430 / 455 / 480 / 505 / 555 / 605 mm	400 / 450 / 500 mm	0° / 3° / 6°	360 / 385 / 410 / 435 / 460 / 485 / 510 mm	7° / 0° - 30° / -6° - 24°	350 - 510 mm
<b>Transit</b>	380 / 405 / 430 / 455 / 480 / 505 / 555 / 605 mm	400 / 450 / 500 mm	0° / 3° / 6°	460 / 485 / 510 mm	7° / 0° - 30° / -6° - 24°	350 - 510 mm

<b>Self propulsion</b>	80° / 0° - 80°	200 - 400 / 330 - 460 / 440 - 510 mm	200 - 280 / 240 / 170 - 300 / 210 - 350 mm	Seat width + 190 mm	305 mm	995 - 1170 mm
<b>Transit</b>	80° / 0° - 80°	200 - 400 / 330 - 460 / 440 - 510 mm	200 - 280 / 240 / 170 - 300 / 210 - 350 mm	Seat width + 160 mm	265 mm	950 - 1050 mm

<b>Self propulsion</b>	850 - 1025 mm	810 - 1040 mm	18°	850 mm	125 / 160 kg*	15,2 kg
<b>Transit</b>	805 - 905 mm	910 - 1040 mm	18°	850 mm	125 / 160 kg*	14,5 kg

<b>Self propulsion</b>	8,3 kg
<b>Transit</b>	11,7 kg

\* Heavy duty version

All product dimensions are in millimetres to apply EN 12183 standards. Our wheelchairs comply with ISO 7176-19 standards and with ISO 10542, and have additional features providing increased levels of occupant security and safety whilst travelling in a motor vehicle. Wheelchair users should however transfer to the vehicle seat and use the vehicle installed restraint system whenever it is feasible. For more pre-sales information about this wheelchair, please consult the [www.invacare.eu.com](http://www.invacare.eu.com) website. It can be found under the section Medium Active Wheelchairs / Action4NG / Sales and Images in the product catalogue. For information about typical user description, intended use, adjustments and transportation, please refer to the user guide.

**Frame colours**



Please notice colors may vary slightly from those displayed above.

Invacare International Sàrl  
Route de Cité-Ouest 2 - 1196 Gland  
Switzerland  
Tel: +41 (0)22 354 60 10  
Fax: +41 (0)22 354 60 11  
[hqeuropa@invacare.com](mailto:hqeuropa@invacare.com)  
[www.invacare.eu.com](http://www.invacare.eu.com)

© 2011 Invacare International Sàrl  
All rights reserved.

All information quoted is believed to be correct at time of print. Invacare reserves the right to alter product specifications without prior consultation.

Action4NG - EU - 06/2011

**Upholstery**



Black Nylon



**Yes, you can.®**